



TECHNICAL SHEET

 **integrale revision**
 **generic revision**

Straight femoral stems product line

Cementless and **Cemented** Version

Ovoid cross section

Same features for both stems




Oblong hole for impaction and threaded hole for extraction, positioned at the top of the stem

Ovoid cross-section for an optimized metaphyseal filling

Rounded edges

Same intra-medullary features for both stems

**More than
20 years
of clinical
experience**

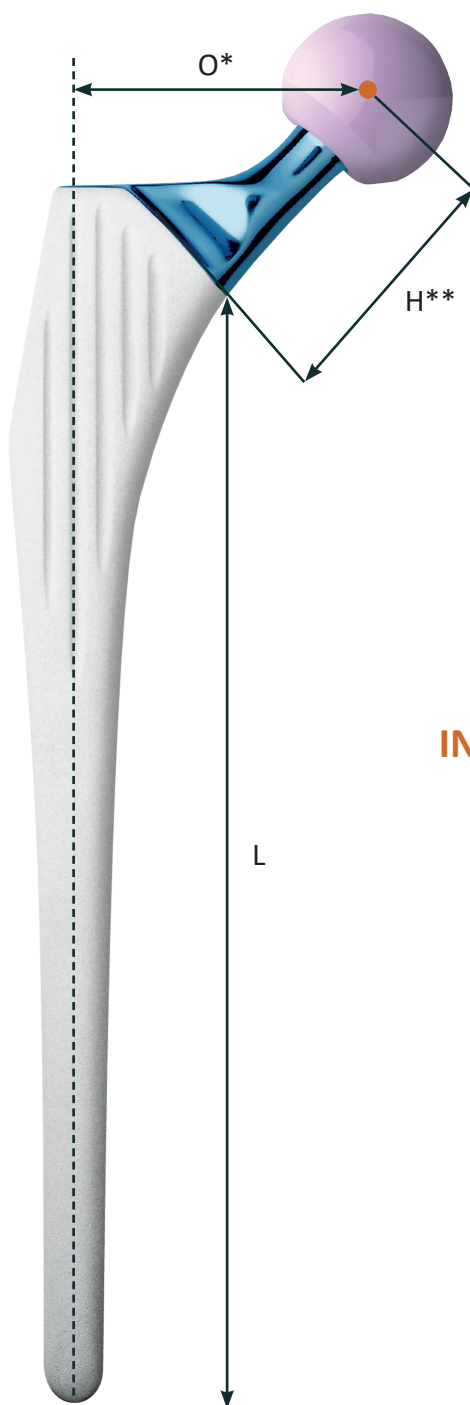

Morse Taper 10/12 (5° 42' 30'')

**Polished/anodised neck
cementless version**

Polished for cemented version.
Neck length: +3.5 mm compared to primary stems of the INTEGRALE range

Restoration of hip biomechanics
Neck-shaft Angle 135°

INTEGRALE /GENERIC Revision

		Offset	Length L
Size 1	O	45.2	182
	H	37.5	
Size 2	O	45.5	187
	H	37.5	
Size 3	O	46.9	192
	H	37.5	
Size 4	O	47.2	197
	H	37.5	
Size 5	O	47.6	202
	H	37.5	
Size 6	O	48.4	207
	H	37.5	
Size 7	O	48.7	212
	H	37.5	

*Distance between diaphysis and head centre with medium neck

**Neck length

integrale revision generic revision

INTEGRALE Revision: Specific Features




Clear visualization of the impaction limit: HA coating

Longitudinal grooves closed on proximal part

HA coating 80 μm


Material: Titanium Alloy (Ti6Al4V) anodized

GENERIC Revision: Specific Features

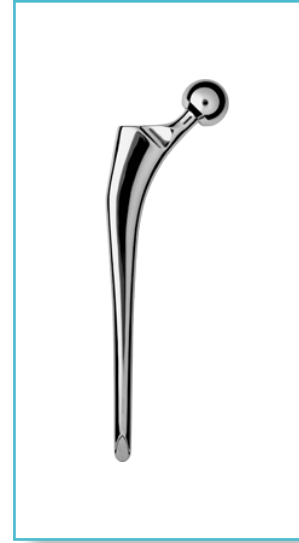

Polished surface on the intra-medullary part of the stem

Material: Stainless Steel (M30NW)

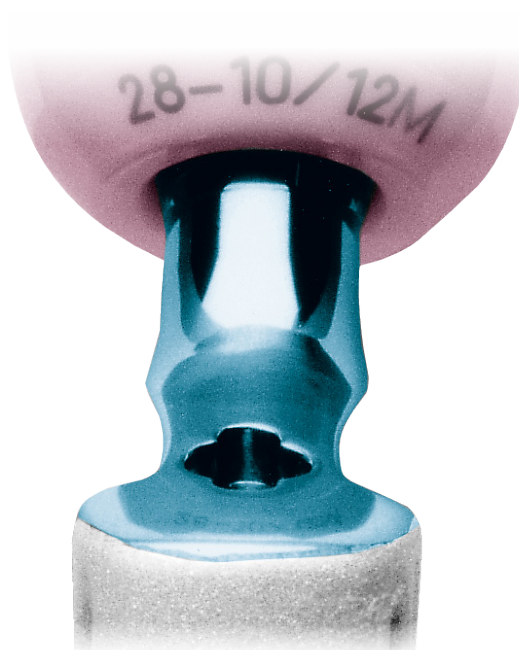

INTEGRALE Revision FEMORAL STEMS

INTEGRALE Revision Femoral Stem - HA Cementless - 10/12 Taper - Size 1	1-01003 01
INTEGRALE Revision Femoral Stem - HA Cementless - 10/12 Taper - Size 2	1-01003 02
INTEGRALE Revision Femoral Stem - HA Cementless - 10/12 Taper - Size 3	1-01003 03
INTEGRALE Revision Femoral Stem - HA Cementless - 10/12 Taper - Size 4	1-01003 04
INTEGRALE Revision Femoral Stem - HA Cementless - 10/12 Taper - Size 5	1-01003 05
INTEGRALE Revision Femoral Stem - HA Cementless - 10/12 Taper - Size 6	1-01003 06
INTEGRALE Revision Femoral Stem - HA Cementless - 10/12 Taper - Size 7	1-01003 07

GENERIC Revision FEMORAL STEMS

GENERIC Revision Femoral Stem - Cemented - 10/12 Taper - Size 1	1-01004 01
GENERIC Revision Femoral Stem - Cemented - 10/12 Taper - Size 2	1-01004 02
GENERIC Revision Femoral Stem - Cemented - 10/12 Taper - Size 3	1-01004 03
GENERIC Revision Femoral Stem - Cemented - 10/12 Taper - Size 4	1-01004 04
GENERIC Revision Femoral Stem - Cemented - 10/12 Taper - Size 5	1-01004 05
GENERIC Revision Femoral Stem - Cemented - 10/12 Taper - Size 6	1-01004 06
GENERIC Revision Femoral Stem - Cemented - 10/12 Taper - Size 7	1-01004 07

FEMORAL heads	Short neck	Medium neck	Long neck	Extra-long neck
M30NW metallic femoral head—10/12 Taper— Ø22,2 mm	1-01011 01 (-2mm)	1-01011 02 (0mm)	1-01011 03 (+2mm)	-
M30NW metallic femoral head—10/12 Taper— Ø28mm	1-01011 11 (-3.5mm)	1-01011 12 (0mm)	1-01011 13 (+3.5mm)	1-01011 14 (+7mm)
M30NW metallic femoral head—10/12 Taper— Ø32 mm*	1-01011 21* (-4mm)	1-01011 22* (0mm)	1-01011 23* (+4mm)	-
BIOLOX Delta ceramic femoral head—10/12 Taper— Ø28 mm	1-01911 01 (-3.5mm)	1-01911 02 (0mm)	1-01911 03 (+3.5mm)	-
BIOLOX Delta ceramic femoral head—10/12 Taper— Ø32 mm	1-01911 11 (-4mm)	1-01911 12 (0mm)	1-01911 13 (+4mm)	-

*optional, available on request only
Compatibility available on IFU

Customer Service – France:

Porte du Grand Lyon,
01700 Neyron – France
Tél. : **+33 (0)4 37 85 19 19**
Fax : +33 (0)4 37 85 19 18
E-mail : amplitude@amplitude-ortho.com

Customer Service – Export:

11, cours Jacques
Offenbach,ZA Mozart 2,
26000 Valence – France
Tél. : **+33 (0)4 75 41 87 41**
Fax : +33 (0)4 75 41 87 42

BIOLOX® is a registered trademark
of the Group CeramTec.

www.amplitude-ortho.com