


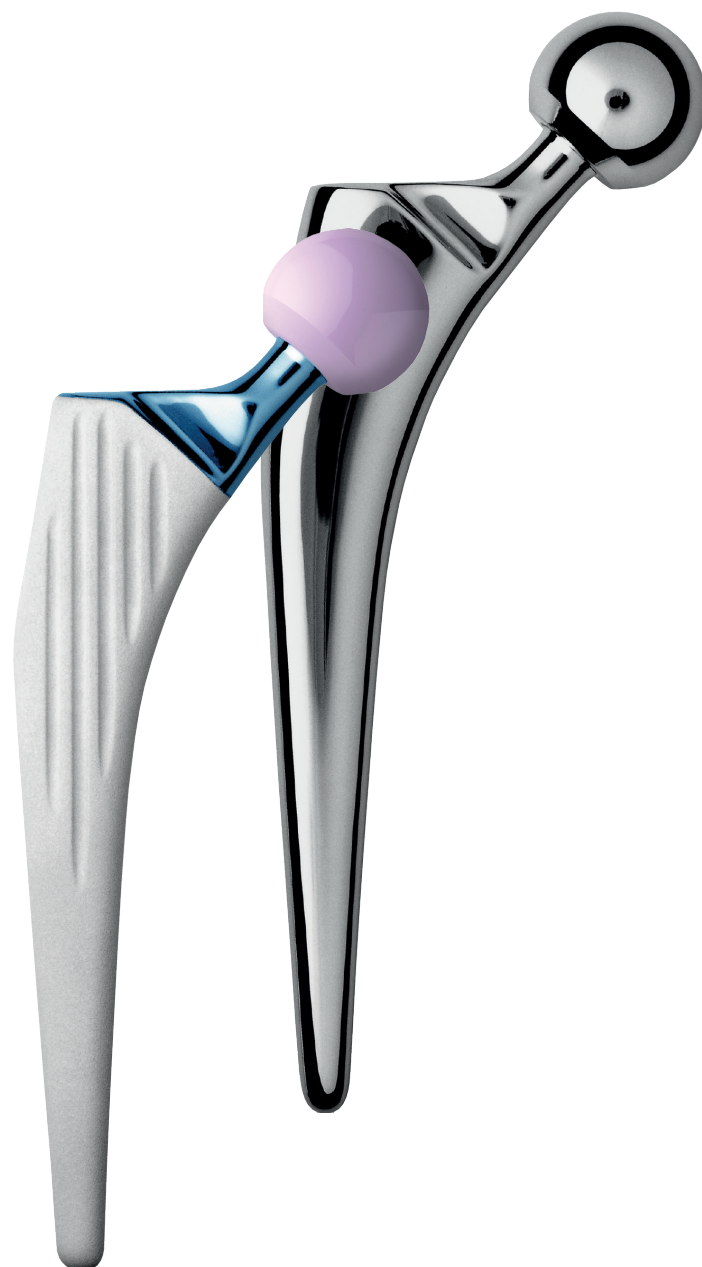

TECHNICAL SHEET

 **integrale**

 **generic**

Straight femoral stems product line

Cementless and **Cemented** Version

Range of 7 sizes with homothetical evolution

Same features for both stems




Oblong hole for impaction and threaded hole for extraction, positioned at the top of the stem

Ovoid cross-section for an optimized metaphyseal filling

Rounded edges

Same intra-medullary features for both stems

More than
20 years
of clinical
experience

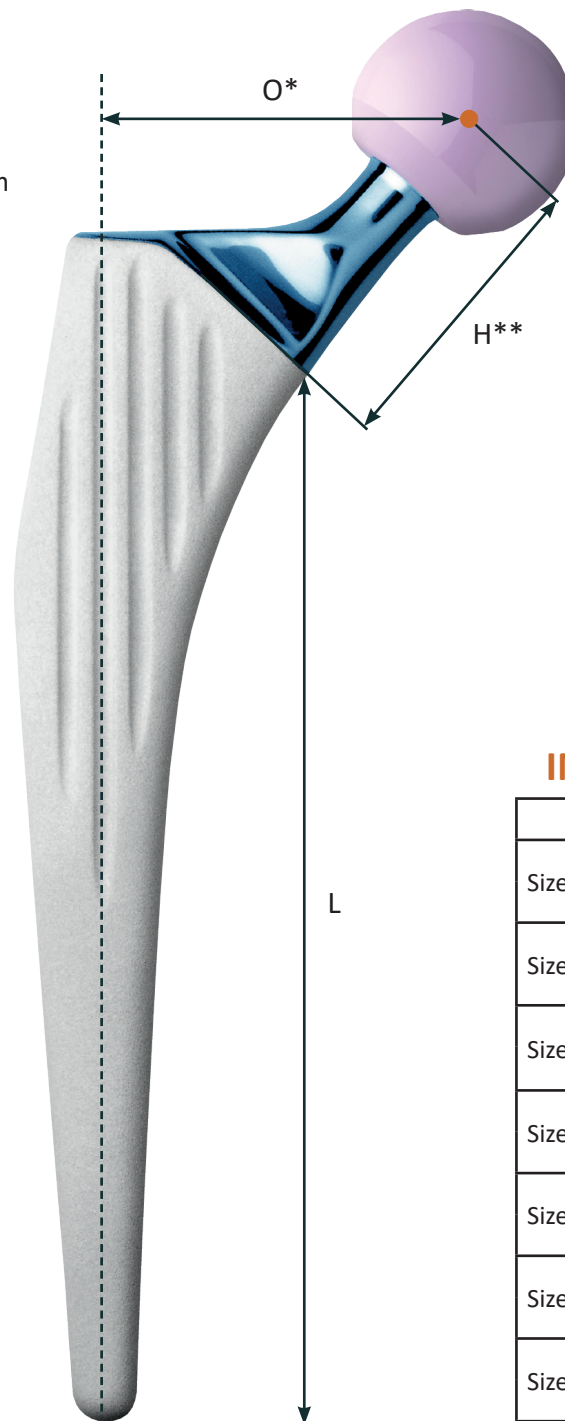

Morse Taper 10/12 (5° 42' 30'')

Polished/anodised neck
cementless version
Polished stem for cemented
version

Restoration of hip
biomechanics
Neck-shaft Angle 137°


INTEGRALE / GENERIC

		Offset	Length L
Size 1	O	41.2	111
	H	34	
Size 2	O	41.7	116
	H	34	
Size 3	O	42.2	122
	H	34	
Size 4	O	42.7	127
	H	34	
Size 5	O	43.2	133
	H	34	
Size 6	O	43.7	138
	H	34	
Size 7	O	44.2	144
	H	34	

*Distance between diaphysis and head centre with medium neck

**Neck length

INTEGRALE: Specific Features

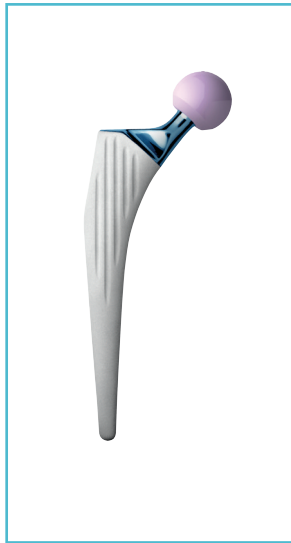

Clear visualization of the impaction limit: HA coating

Longitudinal grooves closed on proximal part

HA coating 80 μ m

Material: Titanium Alloy (Ti6Al4V) anodized

GENERIC: Specific Features



Polished surface of the stem

Material: Stainless Steel (M30NW)

The INTEGRALE system is available on AMPLIVISION:
Assists the surgeon for extra-medullary settings.



INTEGRALE FEMORAL STEMS

INTEGRALE Femoral Stem - HA Cementless - 10/12 Taper - Size 1	1-01001 01
INTEGRALE Femoral Stem - HA Cementless - 10/12 Taper - Size 2	1-01001 02
INTEGRALE Femoral Stem - HA Cementless - 10/12 Taper - Size 3	1-01001 03
INTEGRALE Femoral Stem - HA Cementless - 10/12 Taper - Size 4	1-01001 04
INTEGRALE Femoral Stem - HA Cementless - 10/12 Taper - Size 5	1-01001 05
INTEGRALE Femoral Stem - HA Cementless - 10/12 Taper - Size 6	1-01001 06
INTEGRALE Femoral Stem - HA Cementless - 10/12 Taper - Size 7	1-01001 07

GENERIC FEMORAL STEMS

GENERIC Femoral Stem - Cemented - 10/12 Taper - Size 1	1-01002 01
GENERIC Femoral Stem - Cemented - 10/12 Taper - Size 2	1-01002 02
GENERIC Femoral Stem - Cemented - 10/12 Taper - Size 3	1-01002 03
GENERIC Femoral Stem - Cemented - 10/12 Taper - Size 4	1-01002 04
GENERIC Femoral Stem - Cemented - 10/12 Taper - Size 5	1-01002 05
GENERIC Femoral Stem - Cemented - 10/12 Taper - Size 6	1-01002 06
GENERIC Femoral Stem - Cemented - 10/12 Taper - Size 7	1-01002 07

FEMORAL HEADS	Short neck	Medium neck	Long neck	Extra-long neck
M30NW metallic femoral head—10/12 Taper— Ø22,2 mm	1-01011 01 (-2mm)	1-01011 02 (0mm)	1-01011 03 (+2mm)	-
M30NW metallic femoral head—10/12 Taper— Ø28mm	1-01011 11 (-3.5mm)	1-01011 12 (0mm)	1-01011 13 (+3.5mm)	1-01011 14 (+7mm)
M30NW metallic femoral head—10/12 Taper— Ø32 mm*	1-01011 21* (-4mm)	1-01011 22* (0 mm)	1-01011 23* (+4mm)	-
BIOLOX [®] Delta ceramic femoral head—10/12 Taper— Ø28 mm	1-01911 01 (-3.5 mm)	1-01911 02 (0mm)	1-01911 03 (+3.5mm)	-
BIOLOX [®] Delta ceramic femoral head—10/12 Taper— Ø32 mm	1-01911 11 (-4mm)	1-01911 12 (0 mm)	1-01911 13 (+4mm)	-
BIOLOX [®] Delta ceramic femoral head—10/12 Taper— Ø36 mm	1-01911 21 (-4mm)	1-01911 22 (0mm)	1-01911 23 (+4mm)	-

*optional, available on request only
Compatibility available on IFU

Customer Service – France:

Porte du Grand Lyon,
01700 Neyron – France
Tél. : **+33 (0)4 37 85 19 19**
Fax : +33 (0)4 37 85 19 18
E-mail : amplitude@amplitude-ortho.com

Customer Service – Export:

11, cours Jacques
Offenbach,ZA Mozart 2,
26000 Valence – France
Tél. : **+33 (0)4 75 41 87 41**
Fax : +33 (0)4 75 41 87 42

www.amplitude-ortho.com